


THE TULIP

PRESHAPED EMBRYO TRANSFER CATHETER WITH TULIP SHAPED BULB

The Tulip is a set for difficult cervical passage, consisting of two parts: a guiding catheter (A) and a loading catheter (C).

FEATURES


- preshaped guiding catheter, ideally curved for easy entry when uterus is flexed
- the Tulip tip enables easy atraumatic entry through the cervical canal
- positioner for depth adjustment
- markings on outer tube for easy positioning
- handle with luer lock fitting for a more secure connection to syringes
- reinforced metal support of loading catheter for optimal control
- arrow marking on outer tube for curve indication

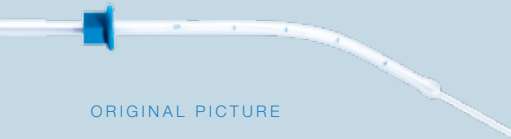

ORIGINAL PICTURE

GYNO INFO THE TULIP

Length	A) 200 mm (\pm 2 mm) C) 266 mm (\pm 2 mm)
Tube OD	A) 2.20 mm (\pm 0.05 mm) C) 1.00 mm (\pm 0.05 mm)
Total volume	C) 0.21 ml
Total volume minus hub	C) 0.069 ml
Tulip bulb	3 mm
Sterile	sterilized by radiation
Packaging	75/shipping box
Quality control	MEA and LAL
Order number	#4000 The Tulip [®]

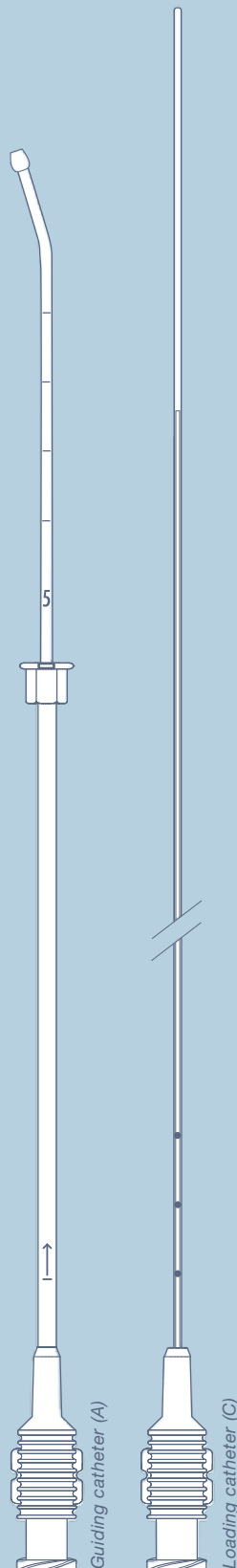

GMP Belgium The Tulip 04/2011



Gynetics Medical Products n.v.
Rembert Dodoensstraat 51 / 3920 Lommel, Belgium
tel. +32 11 64 58 72 / fax +32 11 64 06 77
info@gynetics.be / www.gynetics.be

GYNETICS[®]
Practise better Medicine